

AD 744439

AD

TECHNICAL REPORT

72-35-CE

**FEASIBILITY STUDY TO DETERMINE THE FLUORINE
CONTENT IN QUARPEL-TREATED FABRICS BY ANALYSIS
WITH A 14 MeV NEUTRON ACTIVATION SOURCE**

by

Forrest C. Burns

U. S. Army Materials

and

Mechanics Research Center
Watertown, Mass.

and


Hubertina D. Hogan

and

Gil M. Dias

U. S. Army Natick Laboratories
Natick, Mass.

March 1972

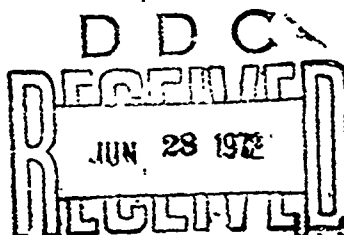

B

Approved for public release;
distribution unlimited.

UNITED STATES ARMY
NATICK LABORATORIES
Natick, Massachusetts 01760



Clothing & Personal Life Support Equipment

Laboratory

TS-181

20

Approved for public release; distribution unlimited.

Citation of trade names in this report does not constitute an official indorsement or approval of the use of such items.

Destroy this report when no longer needed. Do not return it to the originator.

WHITE SECTION <input checked="" type="checkbox"/>		
GOLF SECTION <input type="checkbox"/>		
SPUNDER <input type="checkbox"/>		
SECTION		
DISTRIBUTION/AVAILABILITY CODES		
DIST.	AVAIL.	mm/yr SPECIAL
A		

UNCLASSIFIED

Security Classification

DOCUMENT CONTROL DATA - R & D

(Security classification of title, body of abstract and indexing annotation must be entered when the overall report is classified)

1. ORIGINATING ACTIVITY (Corporate author) U. S. Army Natick Laboratories Natick, MA 01760		2a. REPORT SECURITY CLASSIFICATION UNCLASSIFIED	
		2b. GROUP	
3. REPORT TITLE Feasibility Study to Determine the Fluorine Content in Quarpel-Treated Fabrics by Analysis with a 14 MeV Neutron Activation Source			
4. DESCRIPTIVE NOTES (Type of report and inclusive dates)			
5. AUTHOR(S) (First name, middle initial, last name) Forrest C. Burns Hubertina D. Hogan Gil M. Dias			
6. REPORT DATE March 1972	7a. TOTAL NO OF PAGES 11	7b. NO OF REFS 2	
8a. CONTRACT OR GRANT NO.		8b. ORIGINATOR'S REPORT NUMBER(S) 72-35-CE(TS-181)	
8c. PROJECT NO.			
8d.		9b. OTHER REPORT NO(S) (Any other numbers that may be assigned this report)	
10. DISTRIBUTION STATEMENT Approved for public release; distribution unlimited.			
11. SUPPLEMENTARY NOTES		12. SPONSORING MILITARY ACTIVITY U. S. Army Natick Laboratories Natick, MA 01760	
13. ABSTRACT A study was initiated on an Army-developed water-and oil-repellent fabric known as "Quarpel" containing a fluorochemical component critical to the performance of the finished fabric. To analyze for fluorine content, Quarpel-treated fabrics were subjected to neutron activation. This technique of exposing a fabric to neutrons from a 14 MeV accelerator has not been previously reported. Fluorine content of the treated fabric thus can be measured before and after laundering or dry cleaning and correlated with other physical tests to indicate the durability of the fabric treatment. Details of illustrations in this document may be better studied on microfiche			

DD FORM 1473

REF
OFDD FORM 1473, 1 JAN 64, WHICH IS
FOR ARMY USE

UNCLASSIFIED

Security Classification

UNCLASSIFIED

Security Classification

KEY WORDS	LINK A		LINK B		LINK C	
	ROLE	WT	ROLE	WT	ROLE	WT
Feasibility Studies	8					
Measurement	8,9					
Fluorine	9		7			
Textile Finishes	9		9,4			
Quarrel	9		9,4			
Neutron Activation Analysis	10					
Waterproofing	4					
Oil Resistance	4					
Protective Clothing	4		4			
Physical Tests			8			
Laundering			6			
Drycleaning			6			
ja						

UNCLASSIFIED

Security Classification

Approved for public release;
distribution unlimited.

AD _____

TECHNICAL REPORT

72-35-CE

FEASIBILITY STUDY TO DETERMINE THE FLUORINE CONTENT
IN QUARPEL-TREATED FABRICS BY ANALYSIS WITH A
14 MeV NEUTRON ACTIVATION SOURCE

by

Forrest C. Burns
U. S. Army Materials and Mechanics Research Center
Watertown, Mass.

and

Hubertina D. Hogan and Gil M. Dias
U. S. Army Natick Laboratories

Series: TS-181

March 1972

Clothing and Personal Life Support Equipment Laboratory
U. S. ARMY NATICK LABORATORIES
Natick, Mass.

FOREWORD

This project is a joint effort conducted under the Production Engineering Program, "New Applications for Quarpel." Quarpel-treated fabrics were prepared at the U. S. Army Natick Laboratories prior to analysis for fluorine content using a 14 MeV neutron accelerator at the Army Materials and Mechanics Research Center, Watertown, Mass. The co-authors of this report, Mrs. Hubertina D. Hogan and Mr. Gil Dias of NLABS, performed the laboratory work with the guidance and assistance of Mr. Forrest Burns, the principal author and member of Watertown's Materials Science Division.

The authors wish to acknowledge the cooperation of Dr. Richard N. Macnair, Textile Research and Engineering Division, NLABS, in the preparation of this report.

CONTENTS

	<u>Page</u>
List of Tables	v
List of Figures	v
Abstract	vi
Introduction	1
Technique	2
Experimental	3
Apparatus	3
Procedure	5
Results	5
Evaluation of Untreated Fabric	5
Determination of Positron Emitters in Quarpel-Treated Fabrics	6
Evaluation of Treated, Laundered, Extracted, and Weathered Fabrics	7
Determination of the Homogeneity of Fluorine-Containing Polymer in Commercially Treated Fabrics	8
Discussion	10
Conclusions	11
References	11

LIST OF TABLES

<u>Table</u>		<u>Page</u>
I.	Possible Nuclear Reactions Induced in Quarapel-Treated Fabrics by 14 MeV Neutrons	2
II.	Reproducibility of Fluorine Analysis	8
III.	Percent Fluorine Content of Quarapel-Treated Fabric Samples Containing Two Different Fluorochemicals	8
IV.	Fluorine Content of Fabric Treated with Two Concentrations of Fluorochemical B	8
V.	Fluorine Content of Fabric Treated with Fluorochemical A, Initially and After Laundering and Weathering	9
VI.	Fluorine Content of Commercially Treated Quarapel Fabrics	9

LIST OF FIGURES

<u>Figure</u>		
1.	Neutron Activation Accelerator with Sample Irradiation Holders in Position	4
2.	Sample Irradiation Holder	4
3.	Sample Holders for Counting Irradiated Samples	5
4.	Gamma Ray Spectra Showing the 0.511 MeV Positron Annihilation Peak	6
5.	Exponential Decay Curve Showing Two Different Positron Emitters with Half-lives of 10 and 109 Minutes	7

ABSTRACT

A study was initiated on an Army-developed water-and oil-repellent fabric known as "Quarapel" containing a fluorochemical component critical to the performance of the finished fabric. To analyze for fluorine content, Quarapel-treated fabrics were subjected to neutron activation. This technique of exposing a fabric to neutrons from a 14 MeV accelerator has not been previously reported. Fluorine content of the treated fabric thus can be measured before and after laundering or dry cleaning and correlated with other physical tests to indicate the durability of the fabric treatment.

INTRODUCTION

In 1960, personnel of the U. S. Army Natick Laboratories developed a durable water and oil resistant treatment for textiles to provide better environmental protection for the soldier ^{1/}. This treatment, called Quarpel, is a combination of a quaternary pyridinium salt and a fluorine containing polymer applied in a single bath system. Comparative data showed the superior performance of textiles treated with Quarpel as compared with the standard water repellent finishes then in use. Several variants combining different quaternary compounds and fluorine-containing polymers have since been explored and approved as Quarpel type systems meeting the essential requirements established for the Quarpel label.

A series of tests were devised to evaluate the effectiveness of these various Quarpel treatments. These tests include spray rating, Suter Hydrostatic, dynamic absorption, rain-room simulated rainfall at an intensity of 1"/hour and oil contact angle ratings before and after laundering of the treated textile.

While these tests provide quite complete information on the durability and effectiveness of the Quarpel treatment, a need has been felt for a sensitive analytical method that can reveal the fluorine and nitrogen concentrations on fabrics that have undergone various degradative tests. A non-destructive test technique is preferred since it would permit continued monitorship of a given specimen through a whole sequence of processes and still allow other non-destructive analytical procedures to be made after each process. Chemical analysis and physical test data could thus be compared effectively on the same specimen. This, of course, is not possible with currently used wet chemical analysis techniques.

As a result of discussion with scientists at the U. S. Army Materials and Mechanics Research Center, the use of neutron activation analysis was selected since there are no obvious interfering nuclear reactions (Table I), and a system especially designed to analyze inhomogeneous samples such as Quarpel treated fabrics was available. ^{2/} Fluorine analysis was set as the primary goal since a large amount of nitrogen was already present in some of the fabrics to be analyzed. This was the basis for the cooperative effort, the results of which are reported herein.

TABLE I

Possible Nuclear Reactions Induced in Quarpel
Treated Fabrics by 14 MeV Neutrons

<u>Major Elements Present</u>	<u>Natural Isotope</u>	<u>14 MeV Neutron Reactions</u>	<u>Half-Life</u>
Nitrogen	N-14 N-15	$N^{14}(n,2n)N^{13}$	10 min.
Carbon	C-12 C-13	none none	-- --
Hydrogen	H-1 H-2	none none	-- --
Fluorine	F-19	$F^{19}(n,2n)F^{18}$ $F^{19}(n,p)O^{19}$	109.5 min. 29 sec.
Oxygen	O-16	$O^{16}(n,p)N^{16}$	6.1 sec.

TECHNIQUE

The technique of performing an activation analysis is based on the equation:

$$A = \sigma \phi N(1 - e^{-\lambda t})$$

where

A = number of counts produced by neutron bombardment

σ = nuclear cross section of the reaction in barns

ϕ = neutron flux in neutrons/cm²/sec.

N = number of atoms of the isotope involved in the
nuclear reaction

$(1 - e^{-\lambda t})$ = saturation factor based on the half-life of the
isotope being produced

λ = product decay constant

t = irradiation time (sec)

Actually, by irradiating a standard of known elemental concentration in conjunction with an unknown, σ , ϕ and $(1 - e^{-\lambda t})$ can be cancelled out by saying that the number of counts produced are directly proportional to the molar isotopic concentration of the sample as follows:

$$\frac{A_x}{A_2} = \frac{W_x}{W_2}$$

where

$$W = \frac{\# \text{ of atoms} \times \% \text{ isotopic abundance}}{\text{Molecular wt.} \times 6.023 \times 10^{23} \text{ atoms}}$$

When this equation is modified to apply to this study it becomes:

$$\frac{\%F \text{ of Standard} \times \text{wt. of Standard}}{\text{Counts on Standard}} = \frac{\%F \text{ of Sample} \times \text{wt. of Sample}}{\text{Counts on Sample}}$$

or:

$$\%F \text{ of Sample} = \frac{\%F \text{ of Standard} \times \text{wt. of Standard}}{\text{Counts on Standard}} \times \frac{\text{Counts on Sample}}{\text{wt. of Sample}}$$

EXPERIMENTAL

Apparatus.

The 14 MeV neutron source used in this work was a Kaman Nuclear A-700 Sealed tube unit with an output of approximately 6×10^{10} neutrons/second. The complex irradiation system has been described in detail elsewhere^{2/}, but the accelerator is shown in Figure 1, with the sample irradiation holders in position. These holders are fabricated from polyethylene and have a diameter of 0.35-inches and a length of 1.6-inches (Figure 2). Subsequent to irradiation the samples were taken from the irradiation holders and placed in special counting holders (Figure 3), measuring 1-inch x 7/8-inch x 9/16-inch. The amount of activation was determined by placing the held sample on a 3-inch X 3-inch NaI(Tl) well-crystal* with a Nuclear Data 130A Multichannel Analyzer.

* Thallium-activated sodium iodide crystal.



Figure 1 - Neutron Activation Accelerator with Sample Irradiation Holders in Position

Reproduced from
best available copy.

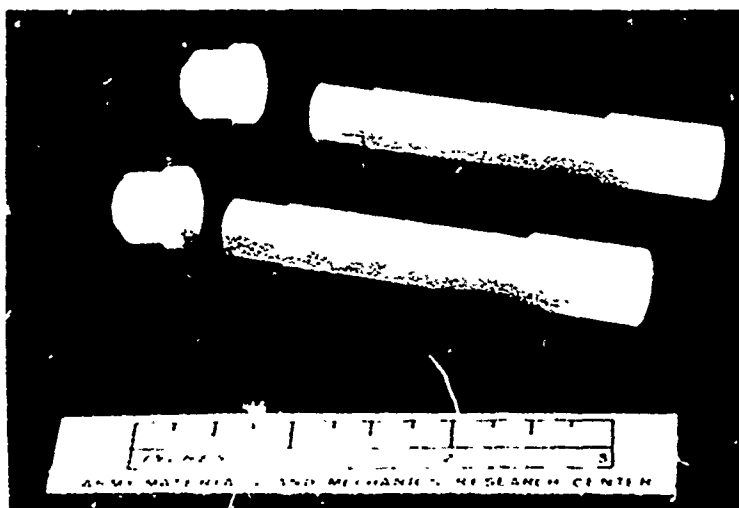


Figure 2 - Sample Irradiation Holders

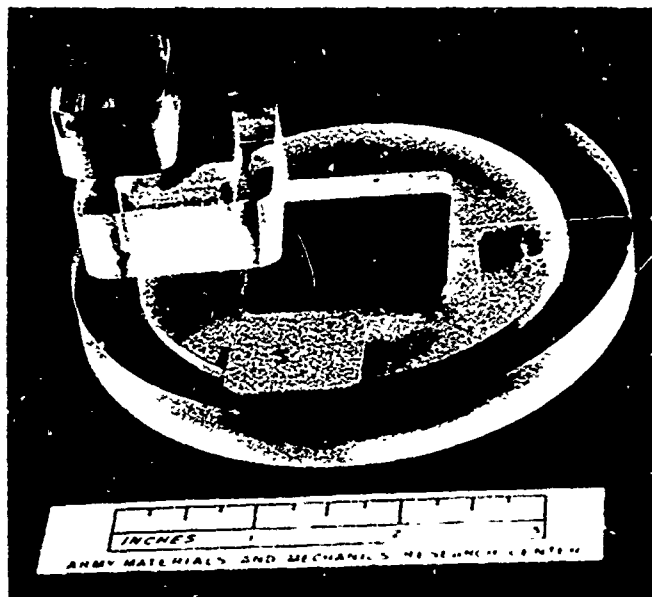


Figure 3 - Sample Holders for Counting Irradiated Samples

Procedure

Samples of the fabric and teflon tape were cut 1-inch x 7/8-inch so that identical geometric shapes were obtained, placed collectively in the irradiation holders, and irradiated simultaneously for a total of 6 minutes. The fabric sample and the teflon standard were then placed in individual counting holders and the amount of irradiation in each determined over a 4 minute period on the NaI(Tl) crystal. The fluorine in the sample was then calculated according to the modified equation given in the Technique Section. The theoretical fluorine content of the teflon tape was used as its actual fluorine content.

RESULTS

Evaluation of Untreated Fabric

Initially untreated fabric samples were irradiated to determine if there were any interfering radioactive nuclides present. The spectra from these samples showed only a positron emitter of 10 minute half-life (Figure 4), indicating the presence of nitrogen, an element known to be in fabrics either as a component thereof (nylon), or in applied treatments (dyestuffs, urea, quaternary nitrogen salts).



Figure 4 - Gamma Ray Spectra Showing the 0.511 MeV Positron Annihilation Peak

Determination of Positron Emitters in Quarrel Treated Fabric

Samples of Quarrel treated fabric were then irradiated and the gamma ray spectra obtained showed only a 0.511 MeV peak on the Nuclear Data 130A Multichannel Analyzer. This indicated only positron emitters to be present. Further study of the exponential decay curve (Figure 5) yielded a composite of two components having half-lives of 10 minutes and 109 minutes; nitrogen ¹³ and fluorine¹⁸ respectively. Since fluorine was the primary interest, the nitrogen ¹³ was allowed to decay out; a period of 2 hours was found to be adequate.

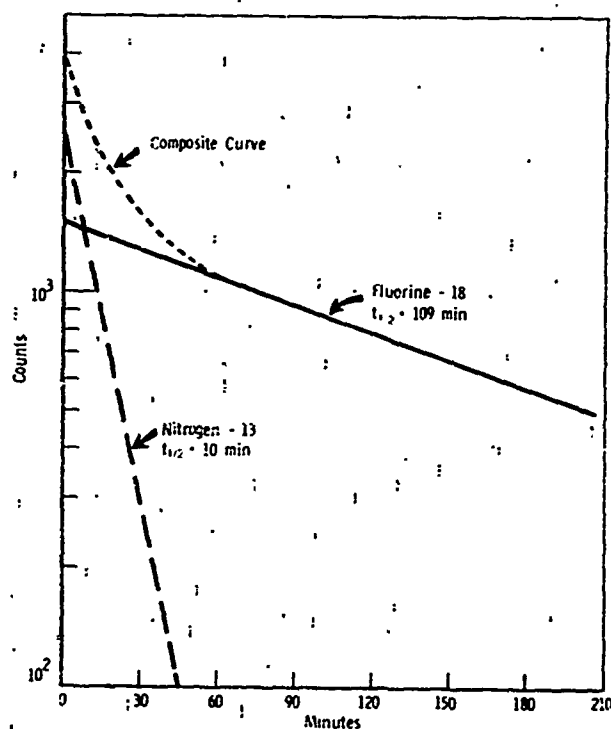


Figure 5 - Exponential Decay Curve Showing Two Different Positron Emitters with Half-lives of 10 and 109 Minutes

Evaluation of Treated, Laundered, Extracted and Weathered Fabrics

Polyester/cotton poplin, 5.7 oz/yd², Cadet Grey 345 fabric was treated in a finishing plant under controlled conditions with two of the accepted fluorine-containing polymers (Product A and Product B). Samples of these fabrics were used to determine:

- A. Reproducibility of the fluorine analysis.
The sample used for this was treated with Product B and the results obtained are listed in Table II.
- B. Comparison of the change in fluorine content of treated fabric containing Products A and B, initially and after laundering and chloroform extraction, (Table III).
- C. The fluorine content of fabric treated with 19.0% and 22.8% fluorochemical B, (Table IV).
- D. The fluorine content of fabric treated with fluorochemical A, initially and after laundering and weathering (Table V).

Determination of the Homogeneity of Fluorine-containing Polymer
in Commercially Treated Fabrics.

Four different fabrics were drawn from the Defense Personnel Support Center: (1) 9 oz/yd² Nyco-sateen; (2) 5.5 oz/yd² cotton oxford; (3) 5.5 oz/yd² cotton warp-nylon fill oxford; and (4) 5.7 oz/yd² cotton/polyester poplin. These materials had been procured by the Army as Quarpel treated fabrics but no further knowledge was received as to their history. The homogeneity was studied by following the fluorine concentration along the width (fill direction) on the first two fabrics and along the length of the last two fabrics (warp direction). The results are shown in Table VI.

TABLE II

Reproducibility of Fluorine Analyses

<u>% Fluorine</u>	<u>Variation</u>
0.64	+0.02
0.59	-0.03
0.62	0.00
0.62	0.00

TABLE III

Percent Fluorine Content of Quarpel Treated
Fabric Samples Containing Two Different
Fluorochemicals

<u>Fluorochemical</u>	<u>Initial</u>	<u>Laundered</u>	<u>CHCl₃ Extracted</u>
A	0.72	0.68	0.44
B	0.70	0.75	0.63

TABLE IV

Fluorine Content of Fabric Treated with
Two Concentrations of Fluorochemical B

	<u>19.0% Fluorochemical B</u>	<u>22.8% Fluorochemical B</u>
	0.64%	0.70%
	0.59%	0.73%
	0.62%	0.69%
Avg.	0.62%	Avg. 0.71%

TABLE V

Fluorine Content of Fabric Treated with Fluoro-
chemical A, Initially and After Laundering and
Weathering

<u>Initial</u>	<u>Laundered</u>	<u>Exposure (Days)</u>		
		18	42	60
0.61%	0.56%	0.63%	0.53%	0.65%

TABLE VI

Fluorine Content * of Commercially Treated Quarpel
Fabrics

<u>Fill Direction</u>	<u>Left End</u>	<u>Center</u>	<u>Right End</u>
	<u>%</u>	<u>%</u>	<u>%</u>
Nyco/sateen, 9 oz/yd ²	0.41	0.55	0.42
Cotton-oxford, 5.5 oz/yd ²	0.78	0.65	0.55
<u>Warp Direction</u>			
Cotton Warp-Nylon Fill oxford, 5.5 oz/yd ²	0.67	0.54	0.52
Cotton/Polyester poplin 5.7 oz/yd ²	0.41	0.56	0.58

* Measurements made at a minimum of 18 inches apart and at
least 3 inches from the selvage.

DISCUSSION

Neutron activation analysis as applied to Quarpel treated fabrics is a new analytical tool. It is a fast, non-destructive method for determining the fluorine content of such fabrics. Results obtained on the AMMRC unique nuclear irradiation transfer apparatus are highly reproducible (Table II). Two of the presently accepted fluorine-containing polymers were used in the Quarpel formulations and compared (Table III). Equal fluorine concentrations were found on both treated samples before and after laundering. However, the sample containing Product A showed greater loss of fluorine when extracted with chloroform than did the Product B sample. Interpretation of this phenomenon can not be made at this time.

This new analytical tool can be used to determine the relative add-on of the polymers used in the Quarpel formulations (Table IV). A slight increase in the percent of Product B in the pad application bath (from 19.0 to 22.8%) produced a 0.09% numerical increase in the percent fluorine content of the treated fabric.

Fluorocarbon loss in laundered and weathered fabrics can be determined by the limited data in Table V. Thus, laundering appears to reduce fluorine content somewhat while weathering causes very little loss on an overall basis.

The inconsistency in the weathering data (Table V) may have been due to inhomogeneity in the sample treatment. Experience has taught that commercially prepared Quarpel treated fabrics had a non-homogeneous finish. However, in general, those fabrics procured under government contracts were able to meet the specification requirements for physical water and oil repellency tests because of the wide tolerance range indicated. Analysis of four stock fabrics drawn from the Defense Personnel Support Center substantiated this hard-earned experience (Table VI) indicating significant non-homogeneity. By effective use of neutron activation analysis these Laboratories will be able to correlate maximum and minimum fluorine contents with relative performance of Quarpel treated fabrics. With such correlations more exact specifications than now possible should be attainable for monitoring procurements, thus assuring the best environmental protection for the soldier that is possible.

CONCLUSIONS

Neutron activation analysis is an effective means of determining fluorine content of Quarpel treated textiles. It is non-destructive and allows correlation of physical test evaluation data with quantitative fluorine data on identical samples. In addition, the decay curve produced (Figure 5) could be resolved with the aid of a least squares computer program allowing the Army to monitor both the fluorochemical and the nitrogen contents of Quarpel treated fabrics and compare them with performance requirements.

REFERENCES

1. "Quarpel Water and Oil Resistant Treatment for Textiles," Textile Series Report No. 111, Headquarters, Quartermaster Research and Engineering Command, Natick, Mass., 1960.
2. Priest, H. F., F. C. Burns, and G. L. Priest, "An Irradiation Transfer and Counting System for Neutron Activation Analysis of Short-lived Components in In-homogeneous Samples," TR 70-21, Army Materials and Mechanics Research Center, Watertown, Mass., 1970.